


# Organic Acid Screening System

## Determination of organic acids in wine and fruit juice

SunData wine system

Method : acids  
Method table : OA0909.tbl  
Sample : new org. acidstd  
Sample number : 369  
Concentration factor : 1.0000  
Standard : c:\sundata\data\S0909336.RPT

**SunChrom**   
...mit uns auf der Sonnenseite der Chromatographie

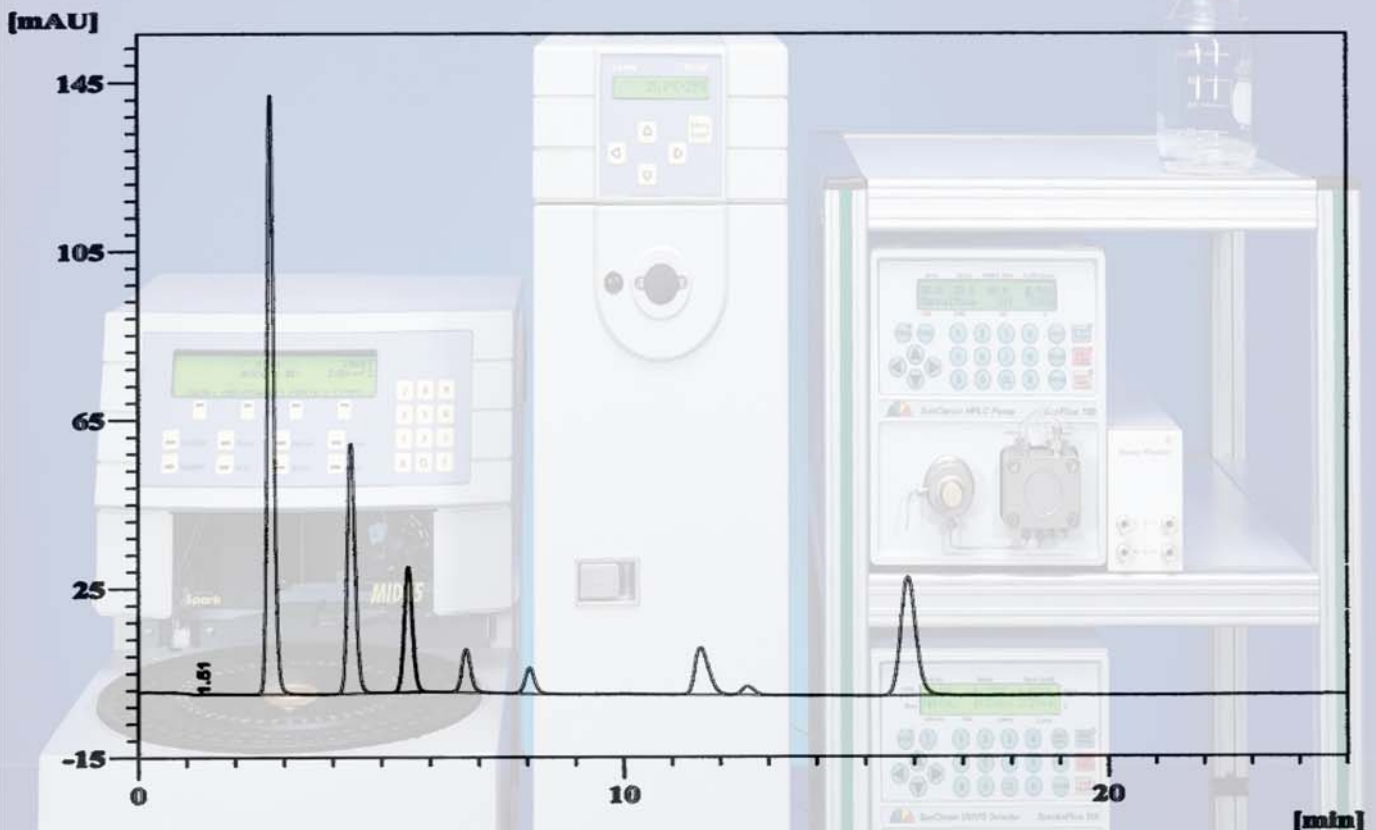

Peak	Retention time	Substance	Peak Height	Amount [g/l]
1	1.51		0.031	
2	2.72	tartaric acid	141.978	4.004
3	4.37	malic acid	59.145	4.345
4	5.54	shikimic acid	29.640	0.055
5	6.72	lactic acid	10.109	2.121
6	8.03	acetic acid	6.001	0.938
7	11.59	citric acid	11.046	1.030
8	12.56	succinic acid	2.040	0.549
9	15.88	fumaric acid	27.977	0.030